

Title : Digestion pattern of pBK549 (Successful)

Date:

03072019

Objective:

To check if the plasmid extracted from the bacteria that contains that plasmid is what we think it is, according to what was sent by Benoit and Agnes from ETH

Method

- Digestion with EcoRv and PvuII enzymes from the general stock
- Using 3ml of bacterial culture from plasmid extraction, using Promega Miniprep Kit

Protocol



Digestion mix

- Final volume of 20ul
- 15,5ul MilliQ
- 2ul plasmid
- 2ul cut smart buffer
- 0,25ul of EcoRv
- 0,25ul of PvuII

Incubation at 37°C for 1 hour, without shaking

Imaging

Results 😊

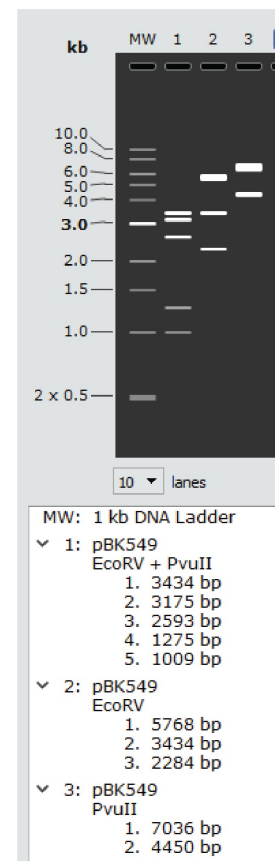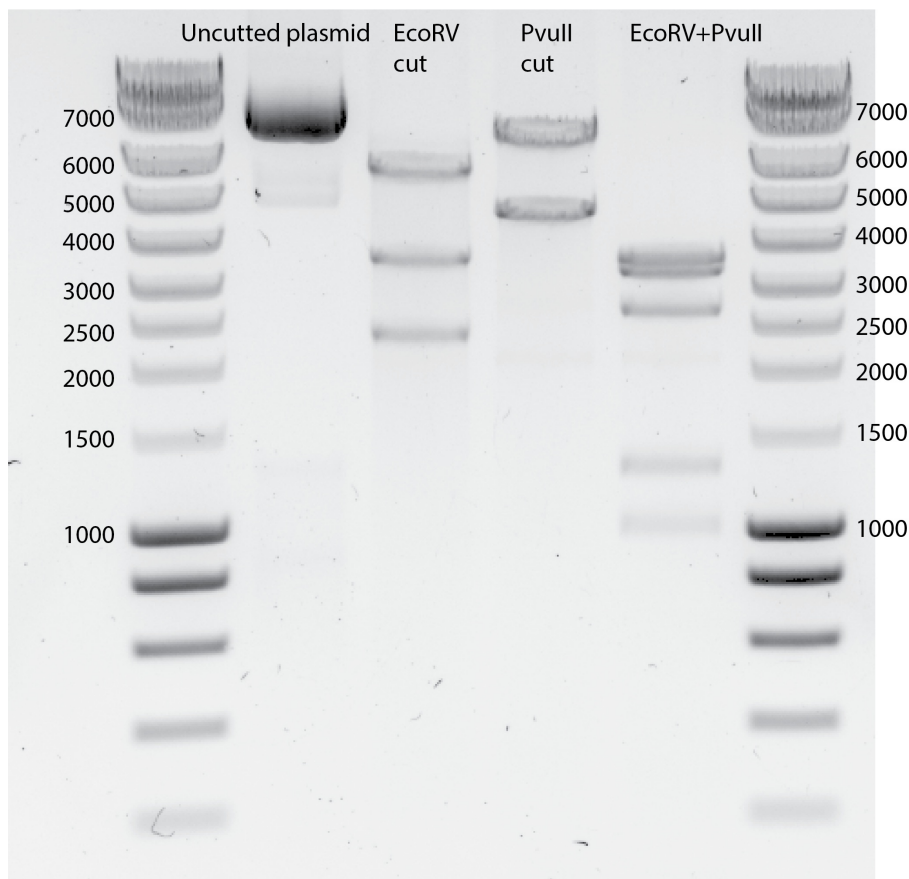

Conclusion